


Stop Patient Sliding – Reduce Cost & Waste

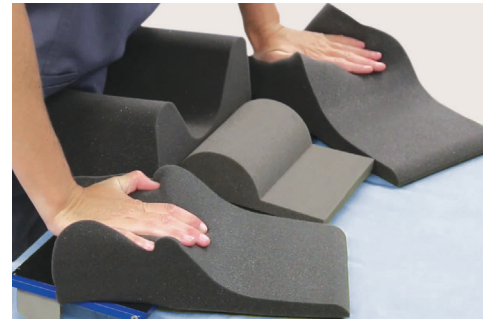

TrenGuard® Trendelenburg Patient Restraint HYBRID

- **No sliding**
Patented bolster keeps patients up to 550 lb (250 kg) from sliding
- **No shoulder braces**
Soft foam lateral stabilizing pillows control body mass shift
- **No chest straps**
Promotes normal ventilation
- **Efficient and compatible**
Accommodates under-patient devices and simplifies patient repositioning

40° Trendelenburg shown with TrenGuard 600 HYBRID Wedge.
TrenGuard Reusable Lateral Stabilizing Pillow ghosted for visualization of bolster.



MADE IN
USA

www.da-surgical.com

info@da-surgical.com

800.261.9953

*compared to viscoelastic static/friction foam pad, shape conforming foam, gel pad, egg crate foam & surgical tape

Pat. www.da-surgical.com/patents

©2024 D. A. Surgical L11147GA7 Rev5

Reduce Waste & Cost-Per-Procedure

TrenGuard® HYBRID Components:




- A** Single-use bolster – prevents patient sliding
- B** Single-use soft head stabilizing pillow – controls side-to-side rocking of the head
- C** Two reusable soft lateral stabilizing pillows – control body mass shift during transition to Trendelenburg

Available in two sizes:

- 1) **TrenGuard 450** for patients up to 450 lb (204 kg)
- 2) **TrenGuard 600** for patients up to 550 lb (250 kg)

Components attach to **D** TrenGuard Patient Support Frame with maximum-strength hook and loop fastener.


Frame secures to O.R. table rails with **E** wide-based rail clamps. TrenGuard is designed to support patients up to 550 lb (250 kg) in extreme Trendelenburg up to 40°.

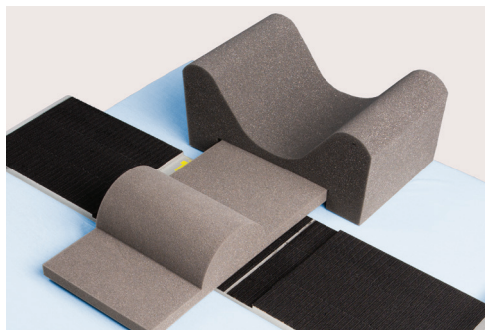

TrenGuard® HYBRID Packs

Each 2-piece single-use pack is vacuum-packed for effective storage (contains bolster and head stabilizing pillow).

Case of 12:
#55201 TrenGuard 450 HYBRID Packs
#56201 TrenGuard 600 HYBRID Packs

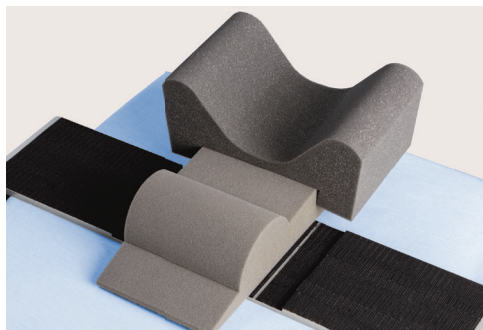

TrenGuard® HYBRID Wedge Packs

For high-BMI patients and those without discernible cervical curve (each pack contains wedge bolster and head stabilizing pillow).

Case of 12:
#55401 TrenGuard 450 HYBRID Wedge Packs
#56401 TrenGuard 600 HYBRID Wedge Packs



TrenGuard® HYBRID Reusable Lateral Stabilizing Pillow

Reduce waste and lower procedure cost. Pillows are washable and reusable up to 500 times.

One pillow:
#55251 TrenGuard 450
#56251 TrenGuard 600

Also available: CLASSIC option with disposable lateral stabilizing pillows, available with standard bolster or wedge; 4-piece single-use packs.

TrenGuard® Patient Support Frame

Item #56100 **D**

Requires two rail clamp (ask about options)

Also available: Bariatric and Dynamic frames (ask for details).



www.da-surgical.com

info@da-surgical.com

800.261.9953

CE USA

Pat. www.da-surgical.com/patents

©2024 D. A. Surgical L11147GA7 Rev5