
Instructions for use:

ShroudGuard™

OR Table Pedestal Protector

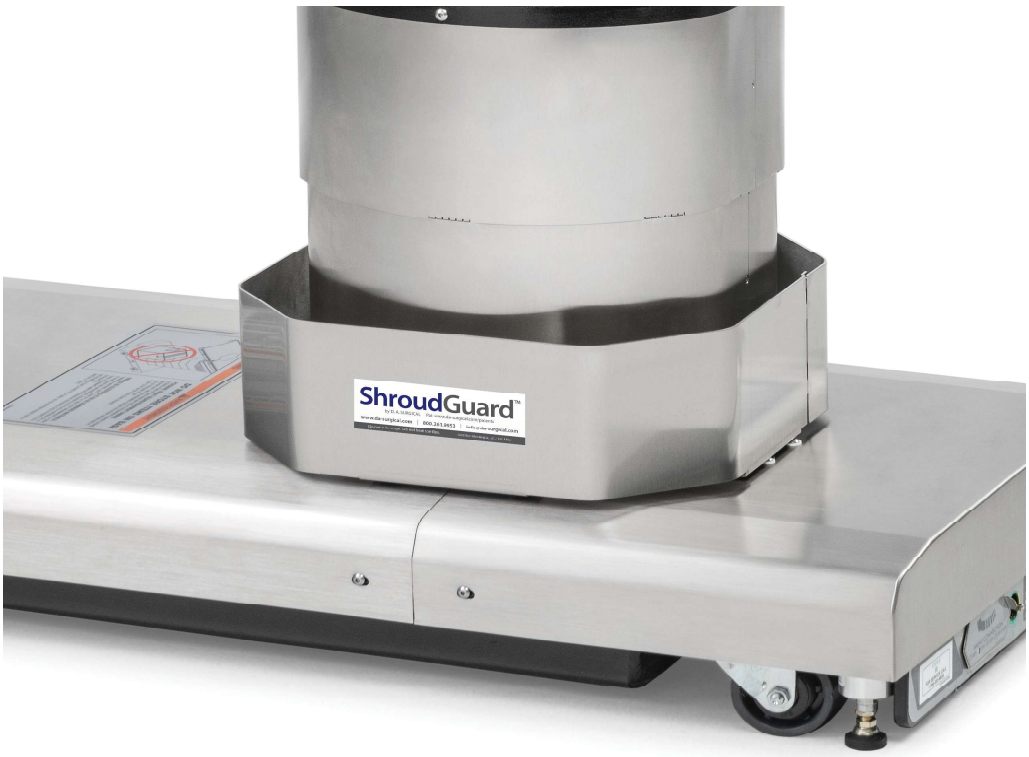

Model#37312 ShroudGuard™ Surround 4085+ shown above.

Rev. 15 6/5/2025

CONTENTS

INTENDED USE..... 2
WARNINGS & CAUTIONS..... 2
INSTRUCTIONS FOR USE..... 3-4
STORAGE & CLEANING..... 5
REPLACEMENT PARTS & SERVICE.... 5
DISCLAIMER..... 5
CE LABEL.....6
EXPLANATION OF SYMBOLS.....6

37110	ShroudGuard™ Surround 7000	OR Table Pedestal Protector for Trumpf® 7000/7000dv
37303	ShroudGuard™ Surround MPB3	OR Table Pedestal Protector for Merivaara Promerix™ and Promerix™ B3
37210	ShroudGuard™ Surround 3085	OR Table Pedestal Protector for STERIS® Amsco 3080, 3085
37220	ShroudGuard™ Surround CMAX	OR Table Pedestal Protector for STERIS® CMAX®
37216	ShroudGuard™ Surround OT1200	OR Table Pedestal Protector for STERIS® OT1200
37310	ShroudGuard™ Surround 3085+	OR Table Pedestal Protector for STERIS® Amsco 3080, 3085
37312	ShroudGuard™ Surround 4085+	OR Table Pedestal Protector for STERIS® 4085, 4095
37314	ShroudGuard™ Surround 5085+	OR Table Pedestal Protector for STERIS® 5085, 5095
37324	ShroudGuard™ Surround OD	OR Table Pedestal Protector for Stryker/Berchtold® Operon D770/760, D820/830, D850/860
37120	ShroudGuard™ Surround Titan	OR Table Pedestal Protector for Trumpf® Titan®
37212	ShroudGuard™ Surround 4085	OR Table Pedestal Protector for STERIS® 4085, 4095
37214	ShroudGuard™ Surround 5085 NL	OR Table Pedestal Protector for STERIS® 5085, 5095 - no latch
37217	ShroudGuard™ Surround 5085	OR Table Pedestal Protector for STERIS® 5085, 5095
37221	ShroudGuard™ Surround CMAX3	OR Table Pedestal Protector for STERIS® CMAX® 3
37222	ShroudGuard™ Surround for CMAX BP	OR Table Pedestal Protector for STERIS® CMAX® 2 with black plastic
37224	ShroudGuard™ Surround OD850	OR Table Pedestal Protector for Berchtold® Operon® D850

D. A. SURGICAL

INTENDED USE

ShroudGuard is a device in the form of a fence/surround intended to protect the pedestal/cladding of an operating table from direct physical damage. It is a multi-piece shield/guard which clips together around the perimeter of the base of the pedestal. This is a reusable device.

INDICATION FOR USE

ShroudGuard should only be used with the appropriate corresponding model of OR table.

WARNINGS and CAUTIONS

⚠ **WARNING - PERSONAL INJURY HAZARD:**


- Read and understand all instructions presented in this manual before attaching/removing ShroudGuard.
- Take care to prevent pinching of clothing or skin when connecting the closure mechanism.

WARNING - INFECTION HAZARD:

⚠
- ALWAYS ensure ShroudGuard is cleaned and dried prior to use.
- DO NOT submerge. DO NOT heat sterilize.

WARNING - EQUIPMENT HAZARD:

- DO NOT use if ShroudGuard is damaged or if shroud cladding on OR table is significantly bent.
- ShroudGuard is not intended to protect table/shroud from C-arm damage.
- ShroudGuard may not protect table/shroud from damage in every situation, care should still be taken when placing extra equipment on table pedestal base and articulating table in extreme positions.
- The table legs should not be tucked under the table past 90° when using the ShroudGuard if your table model is an Operon D770 or D760 (model 37324). Tucking the legs when using the 37324 model for the D770 and D760 tables could lead to damage of the device/table.

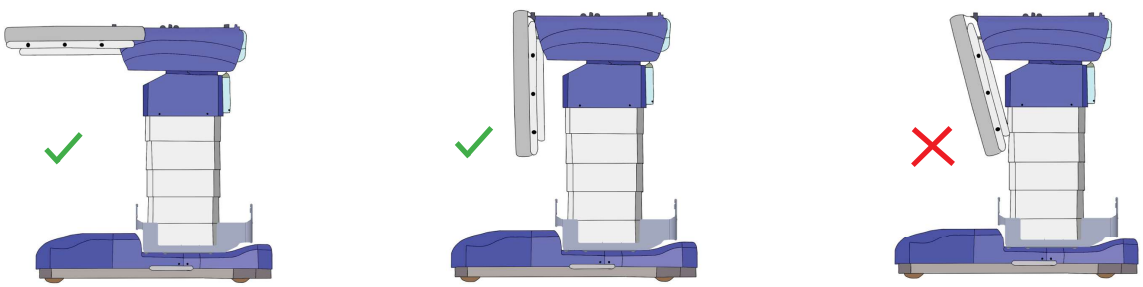

INSTRUCTIONS FOR USE

Step 1 - Verify device and prepare for installation

A. Compare model number and name on UDI label on the ShroudGuard to the listing on the cover of this IFU to ensure you have the correct version for your OR table make and model.

B. Inspect ShroudGuard for damage, warping, bends, dents, or missing rubber bumpers (as applicable) before use.

C. Raise table to full height, with table surface parallel to the floor.


Step 2 - Install

A. Verify proper orientation of ShroudGuard, based on the shape of your table base and/or any labeling or controls on your table shroud cover.

B. Carefully place side 1 around table shroud.

C. Carefully lower side 2 into place, aligning the hooks or tabs with the pin or slots on side 1. Secure ShroudGuard with the applicable closing mechanism as shown below. (Closing mechanism may differ depending on ShroudGuard model number.)


1. Hook Closures: secure as shown in images 1-6 below.
2. Tab Closures: secure as shown in images 7-13 below and on next page.

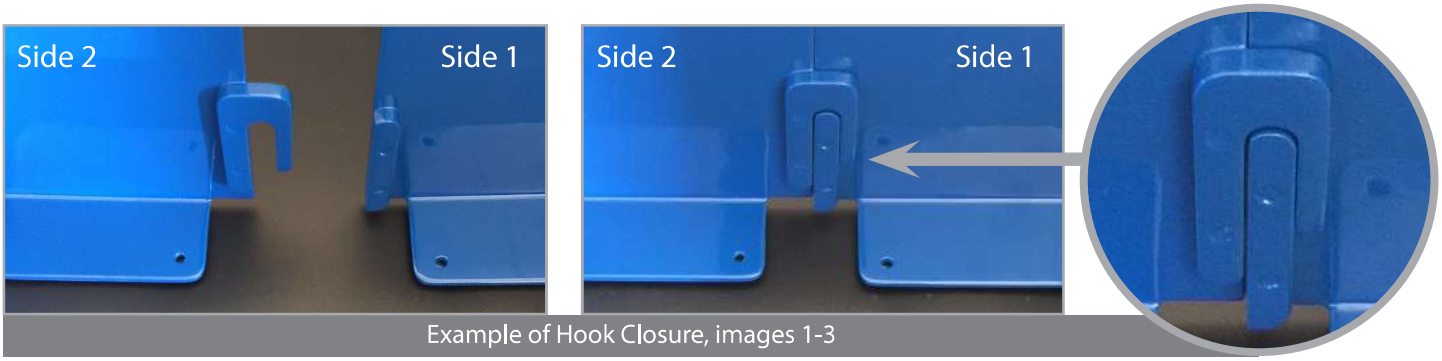

Example of Hook Closure, images 1-3



Example of Hook Closure, images 4-6

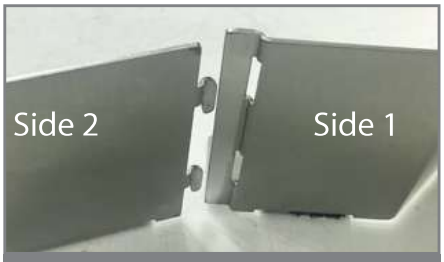

Example of Tab Closure, image 7

INSTRUCTIONS FOR USE (continued)

Step 2 - Install (continued)

3. Tab Closures (continued)
   a. Connect straight tabs into slots of first.




Example of Tab Closure, images 8-10

   b. Swing halves together and connect angled tabs into slots.




Example of Tab Closure, images 11-13

E. If ShroudGuard model uses a Swing Latch Lock, swivel the latch into closed position and confirm it is snug and secure.
(Swing Latch design may differ, depending on ShroudGuard model number.)

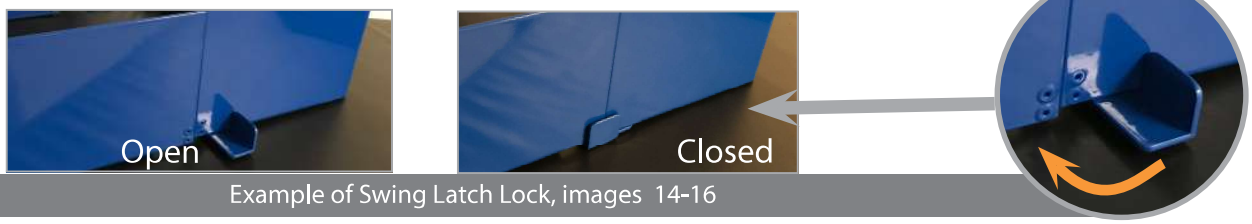

Example of Swing Latch Lock, images 14-16

F. If ShroudGuard model uses a Slide Lock, slide the latch into closed position and confirm it is snug and secure. (Slide Latch design may differ, depending on ShroudGuard model number.)

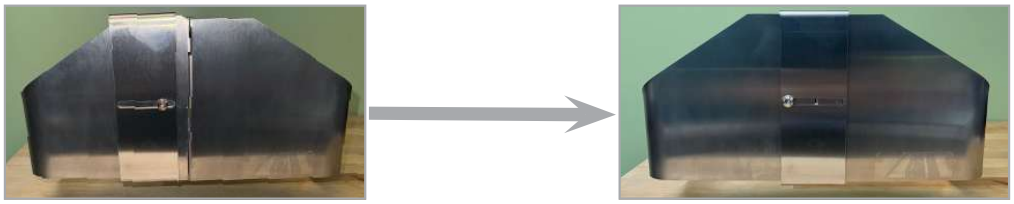

Step 3 - Test

A. Carefully lower the OR table to its lowest setting while checking to make sure that the ShroudGuard is in the proper position and will not make contact with the table shroud.

NOTE: Use caution when articulating OR table in extreme positions to ensure ShroudGuard does not interfere with under-table components or cables.

CLEANING and STORAGE

Follow the hospital's cleaning, hygiene, and disinfection protocol for table accessories. Use hospital-grade cleaners, such as Quaternary Ammonium Compound ("Quats") Germicidal Surface Disinfection/ Deodorizing/Cleaning Wipe or similar to dampen and wipe-down reusable components. Pay particular attention to areas where fluid migration could occur. Allow to dry.

⚠ DO NOT submerge.
DO NOT heat sterilize.

After disinfecting and drying, store ShroudGuard in a secure location to prevent damage or place back into service on appropriate OR table. Take care not to drop device.

REPLACEMENT PARTS and SERVICE

D. A. Surgical warrants that our products will be free from manufacturing defects for a period of 12 months following delivery to the end user.

Useful Life of ShroudGuard OR Table Pedestal Protector is 5 years under normal use.

Always inspect all components before use. Look for signs of warping or cracks/chips in finish which may impact performance.

If unsure of usable condition of ShroudGuard, contact D. A. Surgical immediately at CustomerService@da-surgical.com or 001 (800) 261-9953.

DISCLAIMER

D. A. Surgical assumes no liability arising from misuse or misapplication of ShroudGuard OR Table Pedestal Protector.

It is the sole responsibility of the user and staff to determine applicability of the device for its use, and to study and thoroughly understand these instructions for use and to use ShroudGuard properly.

Contact D. A. Surgical for training options.

D. A. Surgical reserves the right to make changes without notice in design, specifications, and models.


NOTICE

Report any serious incident that has occurred in relation to this device to the manufacturer and the competent authority of the Member State in which the user and/or patient is established

CE LABEL

ShroudGuard
OR Table Pedestal Protector

CE

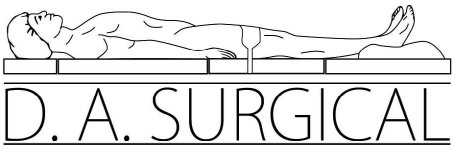

D. A. Surgical
11110 Kinsman Rd Ste 50
Newbury, OH 44065
USA
www.da-surgical.com
customerservice@da-surgical.com

EC REP

Emergo Europe
Westervoortsedijk 60
6827 AT Arnhem
The Netherlands

EXPLANATION OF SYMBOLS

MD MEDICAL DEVICE	
INDICATES A POTENTIAL RISK	CONSULT INSTRUCTIONS FOR USE
LOT BATCH CODE / LOT NUMBER	SN SERIAL NUMBER
MANUFACTURER	DATE OF MANUFACTURE
CE CONFORMITY MARK	EC REP AUTHORIZED REPRESENTATIVE IN EUROPEAN COMMUNITY